


INNOVATION & INTEGRITY SINCE 1984

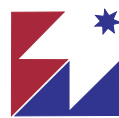

**BOSCO**  
MEDICAL AUSTRALIA

# The **Irwin Re-Breathing Valve** for Adults

*"Feel the Difference"*

The **IRWIN RE-BREATHING VALVE** is the rebreathing valve of choice when it comes to hand ventilating patients during resuscitation, intrahospital transport or physiotherapy.

- **SENSITIVE TO LUNG COMPLIANCE...** Use of the **IRWIN** valve system allows the operator to *feel* changes in a patients lung compliance.
- **GOOD TIDAL VOLUMES DELIVERED AGAINST LOW COMPLIANCE...** With the expiratory flow adjuster fully unscrewed, low resistance to expiratory flow can be achieved.
- **SYSTEM FRIENDLY...** The **IRWIN** valve can be incorporated in most of the Mapleson configurations.
- **USER FRIENDLY...** With a knurled grip on the adjuster and a large button for occluding, the valve can be used in any situation.

#### AREAS OF USE:

ICU  
A&E  
Recovery  
Operating Theatre  
Intrahospital  
Transport

E-mail - [sales@boscomed.com.au](mailto:sales@boscomed.com.au)

**www.boscomed.com.au**

6/70 Andrew Street, Wynnum, Queensland. 4178 FREE CALL: 1800 026 726 Fax: 07 3393 4352

**THE IRWIN VALVE OFFERS THE FOLLOWING ADVANTAGES:**

- Low resistance to expiratory flow.
- Thumb operated button for Inspiratory hold for more uniform tidal volume delivery.
- 2 models - **Tee or Angle** for all situations.
- The flow controlled valve is not pressure controlled by a spring.
- Washable and autoclavable
- Low Resistance



**Part No: IRW01**

**Product Description**

**Part Number**

|  |                  |
|--|------------------|
| Irwin Circuit, straight, silicone - Complete | <b>IRWS001</b>   |
| Irwin Circuit, angled, silicone- Complete    | <b>IRWS002</b>   |
| Irwin Valve - Straight                       | <b>IRW01</b>     |
| Irwin Valve - Angled                         | <b>IRW02</b>     |
| Silicone Tubing - Lightweight                |                  |
| - 30cm                                       | <b>66-10-035</b> |
| - 45cm                                       | <b>66-10-045</b> |
| Rebreathing Bag - 2Lt, (Latex Free)          | <b>825026000</b> |
| Rebreathing Bag - 2Lt, (Silicone)            | <b>66-00-013</b> |
| Circuit Adaptor, 22mm M/M                    | <b>66-12-110</b> |
| Face Mask, Silicone - Size 3                 | <b>3551223</b>   |
| Face Mask, Silicone - Size 4                 | <b>3551224</b>   |



**Part No: IRW02**

**SPECIFICATIONS:**

Weight           99gms  
 Materials:           Stainless & H.D. Plastic  
 Oxygen Flow Rate.   12 to 15lpm.  
 Spring Tension: 3 to 5cm H<sub>2</sub>O @ 15lpm  
 Tapers / connectors: 22mm Male - Tubing end,  
                                   22mm Female - Patient end  
                                   Oxygen, 6mm

**Distributed By:**

**E-mail - sales@boscomed.com.au**

**www.boscomed.com.au**

6/70 Andrew Street, Wynnum, Queensland. 4178 FREE CALL: 1800 026 726 Fax: 07 3393 4352